


# Uricult<sup>®</sup> Vet

In-Office Urine Culture System


*Bring Some Culture to Your Lab*



**Easy to perform • Accurate • Low cost**

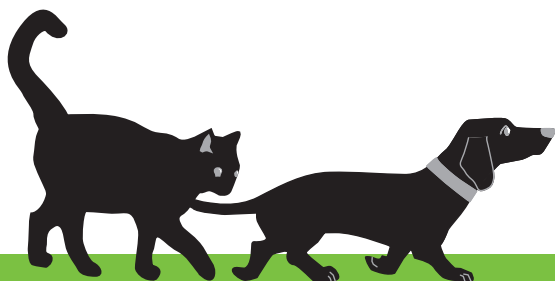

[www.uricultvetusa.com](http://www.uricultvetusa.com)

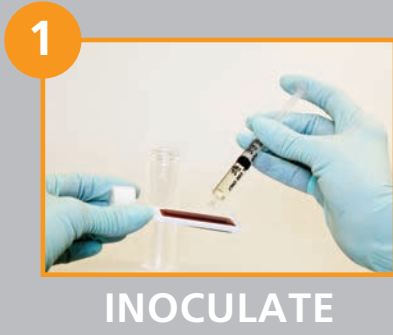

## Uricult® Vet

*Accurate and Easy-to-perform bacterial culture method for the diagnosis of Urinary Tract Infections (UTI)*

The Uricult® Vet dip paddle carries two types of culture media, CLED (green) and EMB (reddish), supporting the growth of the most common uropathogens which cause urinary tract infections in dogs and cats. After inoculation and incubation, the bacteria will form visible colonies on the culture media. Depending on the type of bacteria, the color of the CLED media changes from green to yellow or blue. The result is interpreted by comparing the amount of bacterial growth with the Colony Density Chart provided and using the slide guide or wall chart to make a preliminary identification.

● Uricult® Vet CLED/EMB  
Item # 06254



### Reliable

The test accurately indicates the presence or absence of a urinary tract infection and can provide a preliminary identification of bacterial organisms. The Uricult® Vet dip paddle demonstrates a 99% correlation to the reference laboratory culture method.

### Long shelf-life

Uricult® Vet has a shelf-life of several months at room temperature when protected from drafts and temperature fluctuations. Ideal for any practice size.

### Efficient

Creates a faster and convenient process for veterinary practices and pet owners by identifying a UTI, with accurate results in 16-24 hours. Also, eliminates sending negative samples to reference laboratory while keeping revenue within the practice.

### Great for transportation (if necessary)

Uricult® Vet can be transported without any special precautions to a reference laboratory for further testing and identification (if necessary). Prevents errors in sample handling and maintains the integrity of the specimen during transport.

### Easy to use

Simple three step procedure:

- 1) Inoculate
- 2) Incubate
- 3) Read the results

ORDERING INFORMATION		
PRODUCT	# OF TESTS	ITEM #
Uricult® VET	10	06254
Cultura M Incubator	—	05866



**A PBM Group Company**  
85 Orchard Road • Skillman • New Jersey • 08558  
800.526.2125 • 732.246.3366 • Fax: 732.246.0570 • [www.uricultvetusa.com](http://www.uricultvetusa.com)

Uricult® is a registered trademark of Orion Diagnostica Oy.  
LifeSign is a registered trademark of Princeton BioMeditech Corporation.